

DEBRIDMENT / ROOT CANAL SHAPING

# MicroMega HandFiles

K-Files H-Files K-Reamers




Pre-sterilized instruments  
Blisters of 6 HandFiles

[micromega.coltene.com](http://micromega.coltene.com)

 **COLTENE**

# MicroMega K-Files

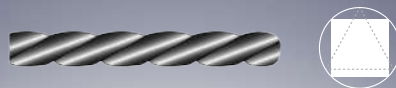

K-Files may be used for manual root canal treatments, as well as for mechanized ones, as valuable supplement to NiTi files to achieve the scouting of the root canal with K-Files before the shaping step, K-Files will help doctors determine the working length, with or without the use of an electronic apex locator, also help remove the primary mechanical constraints through an enlargement action.



Available lengths:

21 mm

25 mm

31 mm



Diameters:

006

to

140



Taper:

2 %


Blister of 6 pre-sterilized instruments, Multiple use

| Lenght | 006      | 008      | 010      | 015      | 020      | 025      | 030      | 035      | 040      | 045      | 050      |
|--------|----------|----------|----------|----------|----------|----------|----------|----------|----------|----------|----------|
| L21    | 60029717 | 60029720 | 60029723 | 60029726 | 60029729 | 60029732 | 60029735 | 60029738 | 60029741 | 60029744 | 60029747 |
| L25    | 60029718 | 60029721 | 60029724 | 60029727 | 60029730 | 60029733 | 60029736 | 60029739 | 60029742 | 60029745 | 60029748 |
| L31    | 60029719 | 60029722 | 60029725 | 60029728 | 60029731 | 60029734 | 60029737 | 60029740 | 60029743 | 60029746 | 60029749 |

| Lenght | 055      | 060      | 070      | 080      | 090      | 100      | 110      | 120      | 130      | 140      |
|--------|----------|----------|----------|----------|----------|----------|----------|----------|----------|----------|
| L21    | 60029750 | 60029753 | 60029756 | 60029759 | 60029762 | 60029765 | 60029768 | 60029771 | 60029774 | 60029777 |
| L25    | 60029751 | 60029754 | 60029757 | 60029760 | 60029763 | 60029766 | 60029769 | 60029772 | 60029775 | 60029778 |
| L31    | 60029752 | 60029755 | 60029758 | 60029761 | 60029764 | 60029767 | 60029770 | 60029773 | 60029776 | 60029779 |

| Assortment | 008/015  | 015/040  |
|------------|----------|----------|
| L25        | 60029921 | 60029841 |
| L31        |          | 60029842 |

# MicroMega H-Files

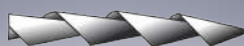

Their steeper-angled cutting edges are best suited for vertical dentin cutting, highly efficient enlargement of the root canal and debris removal from the apex to the coronal part, especially in case of irregular canals walls. Last indication relates to removal of separated instruments.



Available lengths:

21 mm

25 mm

31 mm



Diameters:

006

to

140



Taper:

2 %


Blister of 6 pre-sterilized instruments, Multiple use

| Lenght | 008      | 010      | 015      | 020      | 025      | 030      | 035      | 040      | 045      | 050      | 055      |
|--------|----------|----------|----------|----------|----------|----------|----------|----------|----------|----------|----------|
| L21    | 60029780 | 60029783 | 60029786 | 60029789 | 60029792 | 60029795 | 60029798 | 60029801 | 60029804 | 60029807 | 60029810 |
| L25    | 60029781 | 60029784 | 60029787 | 60029790 | 60029793 | 60029796 | 60029799 | 60029802 | 60029805 | 60029808 | 60029811 |
| L31    | 60029782 | 60029785 | 60029788 | 60029791 | 60029794 | 60029797 | 60029800 | 60029803 | 60029806 | 60029809 | 60029812 |

| Lenght | 060      | 070      | 080      | 090      | 100      | 110      | 120      | 130      | 140      |
|--------|----------|----------|----------|----------|----------|----------|----------|----------|----------|
| L21    | 60029813 | 60029816 | 60029819 | 60029822 | 60029825 | 60029828 | 60029831 | 60029834 | 60029837 |
| L25    | 60029814 | 60029817 | 60029820 | 60029823 | 60029826 | 60029829 | 60029832 | 60029835 | 60029838 |
| L31    | 60029815 | 60029818 | 60029821 | 60029824 | 60029827 | 60029830 | 60029833 | 60029836 | 60029839 |

| Assortment | 008/015  | 015/040  |
|------------|----------|----------|
| L25        | 60029995 | 60029844 |
| L31        |          | 60029845 |

# MicroMega K-Reamers

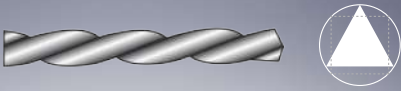

K-Reamers were originally used to find the canals entrances. Today they are mainly used at the end of the root canal preparation as a finisher, to eliminate residual mineral and organic tissues or to place the sealer.



↔

Available lengths:

21 mm

25 mm

31 mm

↔

Diameters:

006

to

140

↻

Taper:

2 %

Blister of 6 pre-sterilized instruments, Multiple use

| Lenght | 006      | 008      | 010      | 015      | 020      | 025      | 030      | 035      | 040      | 045      | 050      |
|--------|----------|----------|----------|----------|----------|----------|----------|----------|----------|----------|----------|
| L21    | 60029654 | 60029657 | 60029660 | 60029663 | 60029666 | 60029669 | 60029672 | 60029675 | 60029678 | 60029681 | 60029684 |
| L25    | 60029655 | 60029658 | 60029661 | 60029664 | 60029667 | 60029670 | 60029673 | 60029676 | 60029679 | 60029682 | 60029685 |
| L31    | 60029656 | 60029659 | 60029662 | 60029665 | 60029668 | 60029671 | 60029674 | 60029677 | 60029680 | 60029683 | 60029686 |

| Lenght | 055      | 060      | 070      | 080      | 090      | 100      | 110      | 120      | 130      | 140      |
|--------|----------|----------|----------|----------|----------|----------|----------|----------|----------|----------|
| L21    | 60029687 | 60029690 | 60029693 | 60029696 | 60029699 | 60029702 | 60029705 | 60029708 | 60029711 | 60029714 |
| L25    | 60029688 | 60029691 | 60029694 | 60029697 | 60029700 | 60029703 | 60029706 | 60029709 | 60029712 | 60029715 |
| L31    | 60029689 | 60029692 | 60029695 | 60029698 | 60029701 | 60029704 | 60029707 | 60029710 | 60029713 | 60029716 |

| Assortment | 015/040  |
|------------|----------|
| L25        | 60029847 |
| L31        | 60029848 |

© COLTENE – www.coltene.com



MicroMega K-Files - MicroMega K-Reamers - MicroMega H-Files  
Advertising | Year of CE-marking: 2020 | Medical device class Is  
according to Directive 93/42/EEC | Notified body: GMED

Medical device for dental care, meant for professional dental use  
only. See product labelling and instructions for use on the website:  
micro-mega.com



Legal manufacturer

**Micro-Mega SA**  
12 rue du Tunnel  
25000 BESANCON / France  
T +33 3 81 54 42 42  
commercial@micro-mega.com

Imprint / Publisher

**Coltène/Whaledent AG**  
Feldwiesenstrasse 20  
CH-9450 Altstätten  
+41 71 757 5300  
info.ch@coltene.com

